


THE SKINCARE TOOLS



LED ROOT RESCUE BRUSH

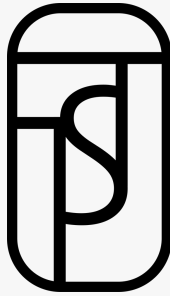


THE SKINCARE TOOLS

At THE SKINCARE TOOLS, we are dedicated to redefining skincare with innovative and sophisticated beauty solutions. By seamlessly blending cutting-edge technology with luxury, we create high-end skincare instruments designed to deliver professional-grade results. Every detail—from the balanced weight of our tools to the premium textures and intuitive designs—has been meticulously crafted to elevate your daily beauty routine into a truly luxurious experience.

The LED Root Rescue Brush is your at-home companion for a revitalised scalp and hair care regimen. Designed to enhance the effectiveness of your favourite hair care serums and oils, it ensures your scalp and strands receive the nourishment they deserve by promoting deeper absorption of nutrients. The integrated LED technology helps improve scalp health and the movement of the comb encourages blood circulation and strengthens hair roots, providing the ultimate foundation for healthier, stronger and more vibrant thick hair.

THE SKINCARE TOOLS

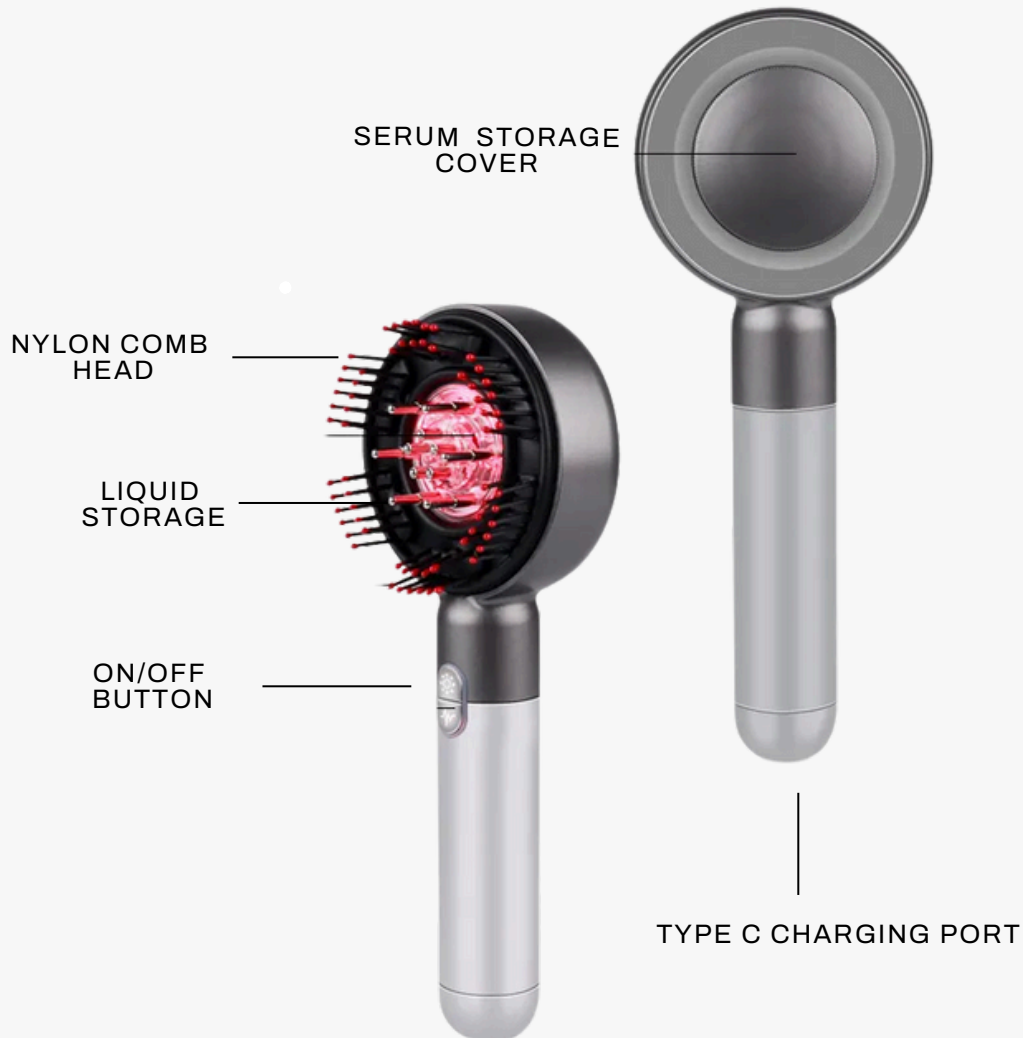

The Serum-Infusing LED Root Rescue Brush is your ultimate at-home companion for enhancing your scalp and hair care regimen. This advanced tool is designed to maximise the effectiveness of your favourite hair care serums and oils, ensuring that both your scalp and every strand receive the nourishment they deserve. By promoting deeper penetration of nutrients, the brush helps to elevate your hair care routine to the next level.

Combining innovative technology with scalp stimulation, the LED Root Rescue Brush boosts blood circulation, strengthens hair roots, and supports healthier hair growth. Its thoughtful design not only improves scalp health but also lays the foundation for vibrant, resilient hair.

FUNCTIONAL INTRODUCTION

Mode Operation

Power On/Off:

- The device features two operation buttons. Hold one of the buttons for 1.5 seconds to turn the device on or off.
- The product cannot be used while charging. After charging, unplug the power adapter before use.

Functional Options:

- Use for 10 minutes each session.
- The two function keys correspond to two modes. To switch modes, press the corresponding mode key while the device is in any working state.
- The machine uses LED lights (red and green) as the gear and charging indicator. In operation, the first gear is green, the second gear is mixed yellow, and the third gear is bright red.
- The indicator light flashes red when charging, and green when full.
- To adjust the gear within the current mode, press the top active mode key.
- Mode gear combination: The default mode combines the massage brush, red LED light, and EMS function combination modes. The red LED light is highlighted by default.
- LED Indicators:
 - The device has two LED colours (red and green) for gear and charging indication.
 - First gear: Green
 - Second gear: Mixed yellow
 - Third gear: Bright red

Serum Addition:

- Open the serum liquid tank cover and add your product as needed. It is recommended to add no more than 3 ml to avoid waste.
- Clean the liquid storage tank after each use to maintain cleanliness.
- Do not add unknown substances or corrosive liquids to prevent damage to the product.

Timing Function:

- Use the device for 10 minutes, 2–3 times per week, or as recommended.
- The device will automatically shut off after 10 minutes of operation.

Product Cleaning:

- To clean the device exterior, wipe it with a damp cloth.
- Rinse the liquid storage bin and liquid outlet hole with clean water after each use to avoid blockages. Dry the product thoroughly and store it in a well-ventilated, dry place.

Storage:

- Store the device in a clean, dry, and ventilated room at a temperature between 5°C and 35°C.
- For long-term storage, pack the product in its original packaging.

FUNCTIONAL INSTRUCTIONS

STEP 1:

TOP OF THE HEAD

- Start from the top of the head.
- Apply a little pressure with the brush on the scalp so the brush teeth have direct contact.
- At the same time, move the brush slowly from the top of the head to the hairline.



STEP 2:

SIDE OF HEAD

- Start from the hair line and fit the brush on the forehead line.
- Move from the middle of the forehead to the ears.



STEP 3:

SIDE OF THE HEAD

- Move in the order of temple-ear-back. Adjust the intensity to desired level.



STEP 4:

BACK OF THE HEAD

- The back of the head is roughly divided into two parts, upper and lower, near the depression. First, stick to the top of the head at the upper part of the depression.
- Move to the lower part of the depression, and start from the back of the neck to the depression of the back of the head and perform the same massage.



PRODUCT FEATURES:

The LED Root Rescue Brush is a LED Red Light Therapy Hair Care Tool that integrates pulse light wave, red light LED and vibration massage functions with a scalp serum applicator. The vibration function gently stimulates the scalp's acupoints, promotes blood circulation, and induces muscle movement, enhancing scalp nutrient absorption and deeply nourishing hair roots.

LED FUNCTION

When the LED light irradiates the scalp, they can stimulate hair follicles, enhance hair firmness and increase hair elasticity, making hair darker and shinier.

VIBRATION FUNCTION

Stimulates the scalp, unclogs pores, activates hair follicles, thereby replenishing hair nutrients and strengthening hair roots. Equipped with a built-in vibration massage device, the vibration mimics the human hand massage, gently touching scalp acupoints, providing more health care for hair and scalp.

SCALP APPLICATOR

The comb's liquid reservoir can be filled with essence, conditioner, etc. The liquid flows out through the comb teeth by pressing, making it easy to comb, nourish hair and scalp, and help strengthen hair roots, rejuvenating hair strands.

NOTE BEFORE FIRST TIME USE

- For the first use, charging is required. Connect the USB cable to charge. After the battery is fully charged, the indicator light will go out. Each charging time is about 1 hour (charging time may vary depending on the remaining power of the rechargeable battery and the usage environment).
- Move slowly on the head and adjust the functions according to your needs for care.
- Avoid direct contact of this product with thicker creams or similar substances.
- Avoid long-term or forceful massage of the same area to prevent skin damage, abrasions, or other issues.
- Avoid looking directly at the LED light for a long time, as it may affect your vision.
- This product is IPX5 waterproof. Do not rinse it under running water or immerse it in water.

BATTERY STATUS

- If no button is pressed, the product will automatically shut down after 10 minutes of use.
- To protect the battery, if the product is not used for more than 2 months, please fully charge it before storing it.

SAFETY PRECAUTIONS

- Avoid direct sunlight, do not expose this product to high temperature or boiling water.
- Please read the instructions carefully and do not use the product continuously in the same part for a long time to prevent causes excessive massage causes skin breakage and scalp pain.
- For hygiene reasons, do not share one product with multiple people. When the product is in working condition, do not leave it alone.
- If you feel any discomfort during use, please stop using it immediately. Keep away from fire and corrosive liquids.
- If the product is damaged, do not continue to use it.
- Do not disassemble, modify, hit or drop this product without authorization. Please place it out of reach of babies, children and pets.
- For your safety, it is strictly forbidden to inject strong acids or liquids harmful to the human body into the product.
- After each use, please pour out the excess liquid in the liquid storage box to avoid solidification. After disassembling the liquid storage box, you can directly rinse it with clean water and then wipe it with a dry cloth.
- When using, you can adjust the vibration and light therapy modes according to your needs for care.
- This product is not for any medical use. If you have the following conditions, it is recommended to consult before use.
- Please do not use any other charging cables to avoid damage of this device.
- Do not use the massage site after surgery, or skin damage, inflammation, congestion and other discomfort symptoms.
- Please avoid dropping or hitting the device to avoid damage.
- **WARNING DO NOT USE IF YOU HAVE THE FOLLOWING:** have medical electronic devices implanted in the body such as pacemakers, suffering from blood circulation disorders or thrombosis, suffering from skin sensitivity or skin infection, heart disease or history, acute diseases, malignant tumors, blood diseases such as bleeding tendency, osteoporosis, soft tissue diseases, stroke and other physical abnormalities.
- **WARNING:** Before using this product, if you do not follow the instructions for operation, The Skincare Tools will not be responsible for any personal injury or economic loss caused.

MAINTENANCE AND SERVICE

- Before cleaning the device, make sure it is turned off and not connected to the power supply.
- After each use, clean the device with a damp cloth and place it in a dry and cool place.
- For personal use only; do not share the device with others; After each use, keep the liquid tank dry and clean.

FAULTS AND REPAIRS

If the device does not work properly, please check the following:

- Make sure the device is fully charged.
- Make sure the "power" button has been pressed and the device is in working mode.
- If there is no response after turning on, it may be that the battery is low, please charge the device and try again.
- This product has an automatic shutdown function after 10 minutes of no button operation.

INCLUDED

- Neo Prene travel bag.
- User Manual.
- Charging Cable.

DISPOSAL OF THE PRODUCT

Proper disposal is essential for environmental safety. When disposing of LED Root Rescue Brush remove the battery first. Ensure the device is turned off before removing the battery. Use appropriate tools to safely open the device without causing damage. Please follow local regulations for the disposal of electronic devices and batteries.

WARRANTY

We appreciate your choice to purchase The Skincare Tools LED Root Rescue Brush. This product comes with a one-year warranty from the date of purchase, covering any technical defects. Please note that the warranty does not cover issues caused by improper use, self repair or damage during transportation.

This device is designed to provide safe and effective at home treatments. Please observe the following guidelines to ensure safe usage and maintain the integrity of the device.